



## ECD004 External Cavity Diode Laser



The MOGLabs External Cavity Diode Laser offers research quality at a teaching-lab price.

The ECD004 is robust and stable. When used with a MOGLabs Diode Laser Controller, mode-hop-free scanning range of up to 40GHz and linewidth below 100 kHz can be achieved using a low-cost *uncoated* 780 nm consumer diode. Diode replacement and re-alignment are easily accomplished by the end-user. Wavelength options extend from 370nm to 1064nm, and powers up to 200mW extra-cavity.

ECD004 is shown above with MOGLabs DLC Diode Laser Controller

### *Features*

- Wide mode-hop free scan range
- Narrow linewidth
- Fast piezo feedback
- Precision alignment controls
- Microwave RF modulation input
- Diode protection circuit and relay
- Low frequency noise
- High feedback bandwidth
- Use with MOGLabs Controller or your electronics

### *Applications*



- Laser cooling and trapping
- Bose-Einstein condensation
- Quantum optics: squeezed light
- Electromagnetic transparency and slow light
- Time and frequency standards
- Laser spectroscopy
- Physics teaching labs

# External Cavity Diode Laser

## Specifications ECD004

Wavelength/frequency	
370nm to 1064nm	Up to 200mW output power, diode dependent
Linewidth	Typically <200kHz, diode dependent
RF modulation	>2.5GHz bandwidth
Coarse tuning range	±5nm for single diode
Optical	
Beam diameter (1/e <sup>2</sup> )	Typically 1mm x 2mm to 1.5mm x 4mm; diode-dependent
Polarisation	Vertical linear 100:1 typical (standard diode)
Thermal	
TEC	±14.5V 3.3A Q = 23W standard
Sensor	NTC 10kΩ standard; AD590, 592 optional
Stability at base	±1mK (controller dependent)
Cooling	Water cooling connections optional (usually not required)
Sweep/scan	
Scan range	Up to 50 GHz, with MOGLabs controller; rate 4Hz to 70Hz
Mode-hop free scan	10 GHz to 40GHz, uncoated diode, with current feed-forward
Piezo	0 – 120V or 0 – 150V, 3 to 5µm
Cavity length	1 – 3cm (5 – 15 GHz FSR) approx.
Electronics	
Protection	Relay, cover interlock connection, reverse diode
Indicator	Laser ON/OFF (LED)
RF input	16MHz – 2.5GHz bias tee (lower cutoff optional)
Connector	MOGLabs DLC Diode Laser Controller (single cable connect)
Dimensions	
Dimensions	105 x90 x90mm (LxWxH), 1kg
Options	

Hermetic seal; water cooling; Faraday isolator; modulation low-frequency cutoff.  
Please contact MOGLabs for further details.

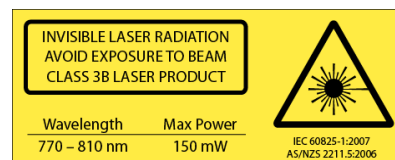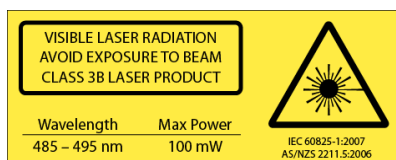

**MOG Laboratories Pty Ltd**  
18 Boase St  
Brunswick VIC 3056  
Australia  
Tel: +61 3 9939 0677  
info@moglabs.com

**MOGLabs USA**  
419 14th St  
Huntingdon, PA 16652  
USA  
Tel: (814) 251-4363  
info@moglabsusa.com

**MOGLabs Europe**  
Goethepark 9  
10627 Berlin  
Germany  
Tel: +49 30 21 960 959  
christoph.p@moglabs.com

© 2015 MOG Laboratories Pty Ltd  
Product specifications and descriptions  
in this document are subject to change  
without notice